



(12) **EUROPEAN PATENT SPECIFICATION**

(45) Date of publication and mention
of the grant of the patent:
23.02.2005 Bulletin 2005/08

(51) Int Cl.7: **G01N 33/551**, G01N 33/552,
G01N 33/547, G01N 33/543,
G01N 33/94, G01N 33/76

(21) Application number: **98303019.8**

(22) Date of filing: **20.04.1998**

(54) **Device and apparatus for the simultaneous detection of multiple analytes**

Vorrichtung und Apparat zur gleichzeitigen Detektion von mehreren Analyten

Dispositif et appareil pour la détection simultanée de plusieurs analytes

(84) Designated Contracting States:
**AT BE CH CY DE DK ES FI FR GR IE IT LI LU MC
NL PT SE**
Designated Extension States:
AL LT LV MK RO SI

(74) Representative: **Perry, Robert Edward**
GILL JENNINGS & EVERY
Broadgate House
7 Eldon Street
London EC2M 7LH (GB)

(30) Priority: **21.04.1997 EP 97302707**

(56) References cited:
EP-A- 0 127 438 **EP-A- 0 768 530**
WO-A-95/16204 **FR-A- 2 693 740**
JP-A- 55 010 546 **US-A- 5 432 099**

(43) Date of publication of application:
28.10.1998 Bulletin 1998/44

(73) Proprietor: **Randox Laboratories Ltd.**
Crumlin, Co. Antrim BT29 4QY (GB)

- **PARSONS R G ET AL: "MULTIANALYTE ASSAY SYSTEM DEVELOPED FOR DRUGS OF ABUSE" CLINICAL CHEMISTRY, vol. 39, no. 9, 1993, pages 1899-1903, XP000676428**
- **HOOK D J ET AL: "SILANIZATION OF RADIO FREQUENCY GLOW DISCHARGE MODIFIED EXPANDED POLY(TETRAFLUOROETHYLENE) USING (AMINOPROPYL)TRIETHOXSILANE" LANGMUIR, vol. 7, no. 1, 1991, pages 142-151, XP000197703**

(72) Inventors:

- **Fitzgerald, Stephen Peter**
Crumlin, Co. Antrim, BT29 4QY (IE)
- **Lamont, John Victor**
Crumlin, Co. Antrim, BT29 4QY (IE)
- **McConnell, Robert Ivan**
Crumlin, Co. Antrim, BT29 4QY (IE)
- **Benchikh, El Ouard**
Crumlin, Co. Antrim, BT29 4QY (IE)

EP 0 874 242 B1

Note: Within nine months from the publication of the mention of the grant of the European patent, any person may give notice to the European Patent Office of opposition to the European patent granted. Notice of opposition shall be filed in a written reasoned statement. It shall not be deemed to have been filed until the opposition fee has been paid. (Art. 99(1) European Patent Convention).

DescriptionField of the Invention

5 **[0001]** This invention relates to a device and apparatus for the simultaneous detection of multiple analytes.

Background of the Invention

10 **[0002]** Traditionally, structurally-diverse analytes have been analysed by means of specific methods, e.g., enzyme immunoassay, high performance liquid chromatography, gas chromatography, enzymatic methods and colorimetric methods. These methods are predominantly one-analyte, one-test methods.

[0003] Automation of analytical methods has generally focused on batch and random access analysers, where multiple analysis on individual test samples is performed using sequential individual test methods. This necessitates the use of multiple packs of individual test kits. In addition, analysis requires employment of several types of equipment, 15 e.g. clinical chemistry analysers, HPLC, GCMS, automated immunoassay instruments or atomic absorption instruments.

[0004] A multi-analyte system should involve a means of providing simultaneous analysis of several analytes in a test sample. This analysis should provide results which identify individual analytes and enable the quantitation of each individual analyte in that test sample. A method of multi-analyte analysis is often claimed but the given criteria are 20 generally not both fulfilled.

[0005] In a multi-analyte system, a typical substrate contains a plurality of individual test reaction sites each possessing a different binding ligand. The test sample contacts each of the reaction zones and thereafter a range of detection techniques is implemented to identify the analyte present. It is important that the detection method used enables quantitation of each individual analyte.

25 **[0006]** In order to produce multi-analyte arrays of spatially-distinct areas of biologically-active ligands on a substrate, the most common approach has been through photolithographic techniques. The substrate is coated with a photolabile linker. In theory, this linker should only become reactive towards a binding ligand following irradiation with light of a suitable wavelength. Spatial resolution is achieved by placing a physical mask (normally manufactured from chrome) on the substrate. The pattern of holes in the mask determines the pattern of binding regions on the substrate.

30 **[0007]** For each biological ligand to be immobilised, the general protocol is: irradiation of first sites, incubation of irradiated substrate with a first ligand to be immobilised, washing to remove loosely bound ligand, blocking unreacted sites activated by the irradiation step, and irradiation of regions where the second biological ligand is to be immobilised, with subsequent steps repeated as for the first ligand. The spatial resolution is dictated by controlling the site of irradiation, either by controlling the site of irradiation by means of a coherent UV light source from a laser or by a number 35 of physical masks and an incoherent light source. This makes the task of immobilising a plurality of biological ligands a time-consuming process. Another disadvantage of the photolithographic approach is the need for expensive physical masks or a laser light source. Further, there is a high degree of non-specific binding.

[0008] For example, the use of arylazides, fluoro-arylazides and benzophenones has been associated with a high degree of non-specific binding. High non-specific binding results in assay background being high, significantly reducing the dynamic range of each multi-analyte assay. The non-specific binding is largely due to passive adsorption of molecules to the non-activated photolabile linker surface through ionic interactions, Van der Waals forces etc.

40 **[0009]** WO-A-9516204 describes a photolithographic approach to reducing the problems associated with high non-specific binding. In this approach, the surface linking molecule was avidin and the photolabile molecule was photobiotin or a derivative thereof. Whilst reduced non-specific binding is claimed, this technique still requires the time-consuming sequences outlined above. Immobilisation of a plurality of 20 separate biological ligands would require a total of 80 steps, assuming the basic requirement of irradiation, binding, blocking and washing steps for each separate ligand to be immobilised.

45 **[0010]** Spatial resolution has also been achieved by passive adsorption. For example, US-A-5432099 discloses binding whereby the ligand to the substrate surface through a combination of ionic interactions, hydrophobic interactions, and Van der Waals forces. Passive adsorption processes are dependent on changes in pH, temperature, ionic strength and on the type of substrate used, making the binding process more difficult to control. The major drawback with this approach is the susceptibility of a proportion of weakly immobilised ligand to be desorbed during the washing step or incubation steps of the biological assay, resulting in poor intra and inter-assay precision.

50 **[0011]** A cross-linker used in many publications has been glutaraldehyde. This linker presents many disadvantages, including the tendency of proteins to cross-link which is likely to alter the function of the protein. A further disadvantage is that the coupling procedure should include a reduction step which is time-consuming and potentially very hazardous, e.g. if sodium cyanoborohydride is used as the reducing agent. Heterobifunctional linkers have been used but in many cases these involve the need for free sulphhydryl groups on the protein to be bound. This necessitates modification of 55

the protein prior to immobilisation.

[0012] In a multi-analyte assay, it is desirable to provide both qualitative and quantitative results. Multi-analyte assays have been available for antibiotics, for example. These are largely based on microbial inhibition assays, where an antibiotic present in the sample inhibits bacterial growth and forms a zone of clearance which is proportional to the concentration of the antibiotic present in the sample. However, this method cannot provide any indication as to the identity of the antibiotic, or an accurate determination of its concentration. Microbial inhibition methods are also very slow, the complete process taking several days.

[0013] Chemical screening methods such as high performance liquid chromatography (HPLC) or gas/liquid chromatography mass spectrometry (GCMS/LCMS) struggle to accommodate the structural diversity/polarity extremes of each antibiotic group, e.g., penicillins, sulphonamides, aminoglycosides, tetracyclines etc. In addition, chromatographic methods necessitate extensive sample preparation in order that the signal-to-noise ratios are such that the required detection limits can be achieved.

[0014] Available multi-analyte devices include the Triage (see Clinical Chemistry 38(9):1678-1684 (1992)) and Advisor (see Clinical Chemistry 39(9):1899-1903 (1993)). These devices are only suitable for qualitative analysis.

[0015] The Triage device is for the simultaneous detection of a panel of seven drugs of abuse in human urine. Each device is only capable of analysing one urine sample. At the end of the procedure, the operator visually examines each of the drug-specific test zones for the presence of a red bar. All steps of the assay protocol must be performed manually by the operator. There is also no hard-copy of the test result available.

[0016] The Advisor device is similar in its application to that of the Triage device. The Advisor device screens for five different classes of drugs of abuse. The device operates using agglutination assay principles, with individual channels for each drug. All steps of the assay protocol are performed by the operator. Negative samples have agglutinated particles, whereas positive drug samples provide an unaggregated pattern of particles.

[0017] The majority of biosensors/microfabricated devices for biological application have employed silicon as the substrate. Others use glass or quartz substrates. Silicon has a very controlled crystallographic structure with well-defined crystal planes. The uniformity of the silicon substrate makes it an ideal choice for the development of a multi-analyte test device.

[0018] However, dark substrates such as silicon give so-called black body effects. In the case of detection by fluorescence, in which a fluorophore is excited using light of a particular wavelength, a dark substrate may absorb the incident excitation light energy, thus diminishing light emission from the fluorophore.

[0019] EP-A-01 27438 discloses the production of a solid state device by photo-activation, using a mask, and bonding a ligand to the thus-activated parts of the substrate surface. WO-A-9516204 also discloses the preparation of a surface-patterned device, by selective irradiation of a photosensitive binding moiety.

Summary of the Invention

[0020] According to the present invention, in a method for the production of a solid state device for performing multi-analyte assays, comprising a substrate and a multiplicity of discrete reaction sites each bearing a ligand covalently bonded to the substrate, wherein the surface of the substrate between the reaction sites is inert with respect to analyte, the method comprises activating, and rendering hydrophobic, the said surface; and applying to the activated, hydrophobic surface an array of ligands, to form the discrete reaction sites, wherein the ligands are applied in water. If desired, the active areas between the reaction sites may be blocked.

[0021] A novel device, obtainable by this method, is a solid state, multi-analyte device which exhibits little or no non-specific binding.

[0022] The ligands may be bound to the substrate via a linker. In particular, it is preferred that the activated surface is reacted successively with an organosilane, a bifunctional linker and the ligand. Bifunctional cross-linkers have been used to provide highly efficient coupling chemistries between organosilanes covalently immobilised on microfabricated silicon or ceramic substrates. Biological ligands can thus be immobilised in multi-analyte arrays.

[0023] This invention removes the need for multiple individual test reagent kits and also several instrument types, thereby facilitating the simultaneous detection of multiple analyses on a single sample. The invention provides a single integrated analyser capable of providing simultaneous detection of a wide range of chemistries. In addition, test reagents may be supplied in a combined format for a particular panel of analytes.

[0024] As part of this invention, a plurality of biological ligands is immobilised in a spatially-defined pattern of spots or lines by means of microfluidic dispensing of the ligand onto a chemically-activated substrate. The biological ligand is covalently attached to the substrate. The coupling efficiency of the biological ligand can be such that the chemical reaction is completed within a few minutes. The immobilisation procedure can ensure that the biological ligand retains its biological activity both in the short term and in the long term.

[0025] This invention also provides an integrated analyser system for the simultaneous detection of a wide range of analytes in a multi-analyte format. The analyser system is designed for maximum end-user convenience, with the

capability of obtaining multi-analyte identification and quantitation from each test sample. The preferred analyser system is a combination of a X-Y translational platform, a sample handling unit, liquid handling/flow control means, a temperature-controlled dark box, a CCD camera, and image-processing software. The platform may be associated with a stepper motor, to achieve a positioned accuracy of, say, 10 μm , for positioning device(s) at each stage of the analytical procedure.

Description of the Drawings

[0026] In the accompanying drawings:

- Fig. 1 shows the formation of a non-uniform substrate surface;
- Fig. 2 shows the chemical activation of groups on the surface of a substrate;
- Figs. 3, 5 and 6 show the covalent immobilisation of a ligand at the surface of a substrate;
- Fig. 4 shows the use of latex particles at the surface of a substrate;
- Figs. 7-11 illustrate devices for the incorporation of chips embodying the invention;
- Figs. 12-14 are schematic views of a system suitable for analysing a device of the invention; and
- Fig. 15 shows calibration curves for analytes assayed by means of the invention.

Description of the Invention

[0027] The substrate that is used in a device of this invention may be, for example, of silicon, quartz, glass or ceramic. Ceramic substrates (aluminium oxide) provide an excellent alternative to silicon substrates, since both fluorescent and chemiluminescent detection techniques can be employed successfully. These findings were unexpected, since the crystallography of ceramic materials would not make them an immediate substrate choice.

[0028] A ceramic substrate may be manufactured to provide a range of grain sizes (1 to 30 μm). The preferred particle size of the ceramic substrate used in this invention is less than 20 μm , preferably less than 10 μm . The reduced particle size imparts much improved surface uniformity which in turn provides enhanced performance of biological assays. Other important features of ceramic substrates include surface topography tolerance, porosity, vacuum-tightness and zero absorption of water.

[0029] The preferred ceramic material consists of 94% alumina (Al_2O_3) with a particle size in the range of 4-8 μm . The material is vacuum-tight, and has a surface topography of 0.6 to 0.8 μm when ground. The surface uniformity can be improved by a polishing process, to yield a surface with variation of 0.4-0.5 μm . A further improvement is achieved by lapping and polishing, to yield a surface with a variability of 0.05-0.1 μm .

[0030] The performance of certain ceramic substrates has been found to be dependent on the characteristics of the grain size. The superior assay performance was found for ceramic materials having a grain size of up to 8 μm , e.g. 4-8 μm . Results for materials of higher grain size were not found to be so satisfactory. Ceramic substrates with grain sizes of approximately 1 μm were evaluated and did not give improved results, with respect to those achieved for the 4-8 μm substrate. This is advantageous, since the material cost is considerably higher (approx. 5-fold) for the very small particle ceramic.

[0031] Suitable silicon substrates are produced with an oxide film of an exact thickness, e.g. a tolerance of $\pm 2\text{nm}$ for a 100 nm oxide film. The oxide film may be 50-500 nm thick, preferably less than 200 nm, more preferably in the region of 100 nm.

[0032] The substrate may be formed as part of a solid-state microfabricated micromachined device, developed for a wide range of panel tests for veterinary and clinical diagnostics applications. Each solid-state test device has an array of reaction regions. Each reaction region is specific for an individual analyte. The reaction region may be in the form of a spot, channel, dimple, pit, well or chamber. The reaction regions are manufactured by immobilising biological molecules onto the substrate.

[0033] Typically, a device of the invention is up to 1 cm^2 in area. The area of each reaction site will usually be less than 1 mm^2 .

[0034] The solid substrate is preferably fabricated to provide an intricate network of ports, chambers, channels, wells, dimples etc. It may also be advantageous to create pillars within the channel or well. Such irregularities can help to achieve maximum surface area interactions between bound biological ligands and test reagents, greatly reducing the incubation times for competitive immunoassays and sandwich immunoassays alike.

[0035] As an alternative or in addition to, say, dimples or channels on the silicon or ceramic substrate, the surface may also be microfabricated to incorporate nanolitre to microlitre mixing chambers/reservoirs/channels. For example, the silicon surface is first oxidised to form an oxide layer. A layer of photoresist is then deposited from which the desired pattern is created. After formation of the pattern on the oxide layer, the photoresist is removed. The silicon is then etched, e.g. using HF, and oxide film removed. Finally, the oxide film is grown uniformly over the entire silicon wafer.

[0036] An illustrative process is shown in Fig. 1. In step (i), a silicon wafer 1 is oxidised to provide an oxide layer 2; in step (ii), a photoresist layer 3 is deposited; in step (iii), the application of light provides a patterned oxide layer; in step (iv), the photoresist is removed; in step (v), the wafer 1 is etched; in step (vi), the oxide film is removed; and in step (vii), a continuous oxide film 2a is reformed.

[0037] The biological ligands are immobilised by covalent bonding. Passive absorption interactions are susceptible to changes in pH, temperature and ionic strength, and may in some instances result in release of weakly-bound ligands during incubation and washing steps, thus contributing to poor assay reproducibility. It is of course desirable that the biological ligand retains maximum activity, after the immobilisation procedure.

[0038] Prior to any chemical activation, the surface of the substrates should be thoroughly cleaned. The first step preferably involves cleaning of the surface by sonication in an alkaline detergent, followed by exhaustive washing with double-deionised water. The substrates are then treated with a chromic acid solution. The chromic acid solution both further cleans the surface and opens surface epoxide groups, as shown in Fig. 2. The epoxide groups may also be opened by other means, e.g. sonication for 1 hour.

[0039] The surface hydroxyl groups thus formed are then available for derivatisation. For example, as shown in Fig. 3, a sequence of reactions comprises the use of an organosilane, then a (hetero)bifunctional cross-linker Z-R-Y, to form a highly reactive intermediate, and finally a functionalised ligand, to give covalent immobilisation.

[0040] In more detail, in the organosilanes of the formula $(RO)_3Si-(CH_2)_n-X$, each R is an alkyl or other hydrocarbyl group such as CH_3 or CH_2CH_3 ; n is an integer, e.g. 1 to 18; and X is a functional group such as epoxy cyclohexyl, NH_2 , CHO, OH, SH, p-chlorobenzyl, m-chlorobenzyl, Br, Cl, $-NH-CH_2-CH_2-NH_2$, 2,3-epoxypropoxy, $-N=C=O$, $-N=C=S$ or p-chlorosulphonylphenyl. These organosilanes may be chosen to provide either a reactive terminal group, capable of forming a covalent bond with a biological molecule, or a less reactive moiety such as NH_2 where further activation with a bifunctional linker is necessary to provide an appropriate end group. Organosilanes possessing terminal electrophilic functional groups do not require activation with a bifunctional cross-linker, since biological ligands can be immobilised covalently through nucleophilic groups on the biological ligand.

[0041] In the case of organosilanes possessing nucleophilic groups, any of a multitude of bifunctional cross-linkers may then be used to provide a very reactive chemical group through which a biological molecule or ligand can be covalently attached. This invention includes the use of bifunctional linkers which can be used in the mass production of chemically-activated substrates and are sufficiently stable to permit long-term storage prior to covalent attachment of the biological molecule or binder ligand. Preferred linkers are inert to normal atmospheric conditions whilst also being sufficiently reactive to form covalent bonds with functional groups of the biological ligand to be immobilised in a very short time period (typically < 10 minutes).

[0042] The bifunctional linker may be, for example, phosgene, thiophosgene, N,N-disuccinimidyl carbonate, xylylenediamine, 1,6-diaminohexane, 1,12-diaminododecane, 1,6-diisocyanatohexane, 1,12-diisocyanatododecane, 1,4-phenylenedithioisocyanate, cyanuric chloride, terephthaldehyde, p-nitrobenzoyl chloride, sulfanilic acid, 2-fluoromethylpyridinium p-toluenesulfonate, 3-aminophenylboronic acid, p-bromophenylboronic acid, diethyl pyrocarbonate, ethyl chloroformate, p-bromoaniline, p-bromophenyl hydrazide, p-bromobenzaldehyde, the 1,2-ethylene glycol of p-bromobenzaldehyde, N,N'-carbonyldiimidazole, terephthaloyl chloride, epichlorohydrin, 1,4-diodobenzene, 1,4-dibromobenzene or a N-hydroxysuccinimide derivative, e.g. of p-aminobenzoic acid, p-bromobenzoic acid, p-bromophenylacetic acid, p-bromoethylbenzoic acid, p-bromomethylbenzoic acid, p-formylbenzoic acid, p-hydroxymethylbenzoic acid, 1,2-ethylene glycol of p-formylbenzoic acid, p-bromophenylpropionic acid or p-hydroxyphenylpropionic acid. A photolabile cross-linker may be used to react with an organosilane having a nucleophilic or electrophilic terminal group. The cross-linker is, for example, the N-hydroxysuccinimide of p-azidobenzoic acid or p-aminobenzophenone.

[0043] Instead of, or in addition to, the use of bifunctional linkers, it is also advantageous to covalently immobilise a layer of latex particles. The diameter of the latex particles is preferably less than 500 nm, and more preferably less than 150 nm. The latex particles may have range of functional groups, e.g. $-CH_2Cl$, $-CHO$, p-chlorophenyl, p-chlorostyryl, $-N=NH$, $-NH-NH_2$ or $-NH_2$. The latex particles may be incubated at a concentration of approximately 0.5 to 1 % w/v with substrates modified with the appropriate organosilane, with or without the presence of a bifunctional linker as described above.

[0044] Fig. 4 shows two reaction schemes for the immobilisation of latex. Either may be followed by activation of the latex with a second linker, and immobilisation of a biological molecule, or direct covalent immobilisation of, say, an antibody.

[0045] An alternative to the use of polystyrene latex particles is the covalent immobilisation of biological ligands to polyethylene glycol derivatives already anchored on a silanated substrate. For example, PEG derivatives with two electrophilic groups such as epoxy or carbonylimidazole are reacted with a silane having a terminal $-NH_2$ group, such as APTES, on a substrate of choice. A suitable reaction sequence is shown in Fig. 5.

[0046] It may also be advantageous to covalently immobilise the biological ligand directly to the silane layer, thus avoiding the need for activation of the silane with polymeric materials or bifunctional linkers. The organosilanes should be receptive to nucleophilic attack by a chemical group (e.g. NH_2 , SH, OH or $NH=NH_2$) on the biological ligand. Orga-

nosilanes that are suitable for direct biological ligand attachment may possess halide, epoxy, isocyanato, aldehyde or tosylate functional groups. Such a reaction is shown in Fig. 6, wherein E is an electrophilic group on the organosilane. Examples of E are Br, Cl, $-O-CH_2-CH=CH_2$, $-NCO$, $-CHO$ and p-chlorosulphonylphenyl.

5 **[0047]** Organosilanes possessing electrophilic groups, e.g., glycidoxy, also have the advantage of being less susceptible to polymerisation during the silanation procedure due to the absence of nucleophilic groups available for attacking the methoxy or ethoxy function of the organosilane. Therefore, the substrate surface should not contain polymerised organosilane.

10 **[0048]** The chemistry of the surfaces provides a means of achieving spatial resolution by virtue of the rapid kinetics of the formation of covalent bonds between the surface chemical functional group and a suitable chemical group present in a sterically favourable position on the biological molecule to be immobilised. The biological molecule is preferably presented to the surface of the substrate by a microfluidic dispenser in the form of an individual droplet or series of droplets which form a line. The volume dispensed is of the order of 1 to 100 nl, preferably less than 50 nl, e.g. closer to 10 nl. The rapidity of the formation of the covalent bonds is such that covalent immobilisation is achieved in minutes, before the dispensed droplet or line evaporates on the surface. The positional accuracy to which the droplet or line of liquid is delivered should have a tolerance of $\pm 20 \mu\text{m}$.

15 **[0049]** The ligands are applied in water, and the surface is hydrophobic, to prevent any lateral diffusion of the dispensed droplet or line. This property contributes to excellent spot quality and reproducibility and enables a greater number of spots of biological molecules to be covalently immobilised per unit surface area.

20 **[0050]** The present invention overcomes the problems associated with conventional photolithography, by enabling the formation of spatially distinct spots of biological ligands, with no requirement for UV light or physical masks. As indicated above, spatial resolution may be achieved by microdispensing techniques. An important factor is the rapid kinetics of the covalent coupling reaction, to ensure highly efficient coupling of the biological ligand in a spatially distinct region, eliminating the lateral digression of the immobilised biological ligand.

25 **[0051]** Unreacted chemical moieties on the substrate may then be blocked, e.g. using blocking molecules known to those skilled in the art. Suitable such molecules include proteins such as casein, bovine serum albumin, lactalbumin etc. or low molecular weight blockers such as glycine, glutamine etc.

30 **[0052]** Photolabile linkers can also be used. For example, the organosilane on the surface of the substrate is reacted in the dark with a photolabile linker (e.g. benzophenone, arylazides etc.) The surface is then spotted with the biological ligands as desired, and covalent attachment is achieved following a short period of irradiation with UV light or a longer period with visible light. The remaining regions of the substrate surface are blocked, using blockers similar to those described molecules as above, in the presence of UV or visible light.

35 **[0053]** The substrate-immobilised biological molecules may be stabilised, e.g. by incubation in a sugar solution (e.g. trehalose) for a short time (1 hour), followed by drying at 37°C for 16 hours. The stabilised substrates may then be sealed in foil pouches with desiccant and stored. The immobilised biological molecules are stable for more than 6 or 12 months, e.g. up to and beyond 2 years when stored at $+2$ to $+8^\circ\text{C}$.

40 **[0054]** The devices may be placed in a range of different carriers which incorporate features which control the efficiency of mixing of test reagents. The flow of liquid test reagents may be achieved by capillary attraction, centrifugal force, vacuum force or electroosmotic flow. The use of electroosmotic flow may avoid the need for valves, so that no moving mechanical parts are used.

45 **[0055]** Closed channels may be formed by bonding a glass plate to the microfabricated surface. The biological molecules are covalently immobilised on the surface prior to bonding the glass plate. Many bonding procedures, e.g. anodic bonding, involve elevated temperatures that may destroy a biological molecule. Therefore, bonding techniques should be non-denaturing, to immobilised biological molecules. One suitable method is indirect bonding, e.g. where the wafer is bonded to a glass plate by a suitable glue, e.g. epoxy glue.

[0056] The devices may then be placed in a suitable carrier. Various such carriers are illustrated in Figs. 7 and 8. By way of example, the dimensions of the device shown in Fig. 8 are 48.62 mm x 48.62 mm, including wells having an internal diameter of 10 mm and an external diameter of 12.82 mm. The centre-to-centre spacing of the wells is 15.36 mm.

50 **[0057]** The devices may incorporate features to enhance mixing of test reagents, samples etc. This is illustrated in Figure 9, where the device includes a reagent addition site 11, reagent channels 12, and test reaction sites 13.

[0058] In Figure 10, the device includes reagent reservoirs 21, a delivery manifold 22, reagent delivery channel test sites 23, and a waste reservoir 24. Fig. 11 shows a 4-channel test structure with similar parts.

55 **[0059]** This invention provides a completely integrated system for the simultaneous, quantitative detection of analytes of a wide range of molecular weights, structural diversity and polarity. Analyte panels are available as appropriate for clinical/veterinary diagnosis or drug screening.

[0060] Depending on the analytes, binding ligands are chosen accordingly. This is within the skill and knowledge of those in the art. Suitable analytes include:

antibiotics, e.g., tetracyclines, sulphonamides, ionophores, aminoglycosides, penicillins or fluoroquinolones;
 hormones, e.g. Luteinising Hormone (LH), Prolactin (PL), Follicle Stimulating Hormone (FSH) or Thyroid Stimulating Hormone (TSH);
 markers of cardiac damage, e.g., myoglobin, carbonic anhydrase, troponin I, glycogen phosphorylase BB, CK-MB,
 5 fatty acid binding protein or troponin T;
 markers of infectious disease;
 allergy markers;
 drugs of abuse;
 enzymes;
 10 viruses;
 nucleotides; and
 peptides.

[0061] For example, one panel is for the detection of sulphonamide antibiotics. This invention provides a method for
 15 the simultaneous quantitative identification of, say, up to 20 individual sulphonamides. Other examples include cardiac, fertility and infectious disease panels.

[0062] The invention makes it possible to simultaneously detect up to, say, 20 analytes which may have no structural
 similarity. Sample matrices that may be tested include serum, plasma, urine, bile, faeces, tissue, water and feed. The
 20 volume of sample required is very low, typically < 1.5 µl/analyte. The test reagents, e.g. enzyme-labelled antibodies, enzyme-labelled haptens, fluorescently-labelled antibodies or fluorescently-labelled haptens, may be all contained in a single reagent reservoir, dramatically reducing the liquid-handling requirements.

[0063] In sandwich assays, e.g. of Luteinising Hormone, Follicle-Stimulating Hormone, prolactin, Thyroid-Stimulating
 Hormone etc, the sample is added along with an assay buffer and incubated for a short period which is typically less
 25 than 30, and preferably less than 10, minutes. Following a wash step, the cocktail of labelled detecting antibodies is added and incubated for a further period of time. This period is again typically less than 30, and preferably less than 10, minutes. The device is then washed, to remove any unbound label, and the signal quantified.

[0064] It may be advantageous for certain assays to incorporate a facility within the microfabricated device to remove
 potential interferents such as rheumatoid factor interference. Removal of rheumatoid factor may be achieved by con-
 30 tacting the test sample to an area of immobilised immunoglobulin, for example prior to the test sample contacting the reaction region.

[0065] A further example is HAMA (Human Anti-Mouse Antibodies) interference; these antibodies can cause severe
 problems in the performance of assays utilising monoclonal mouse antibodies. The traditional solution is to include
 expensive additives in test reagents to counteract the problem. In this invention, there is the advantage of removing
 35 the HAMA interference by contacting the sample with regions on the microfabricated device, to remove these antibodies, prior to the sample reaching the reaction region.

[0066] More generally, ligands may be provided over part of the device, that bind contaminants. This is especially
 valuable where defined spreading is allowed on the surface of the device, e.g. in channels. The capability of removing
 components that interfere enhances the accuracy of the results generated.

[0067] The detection labels may also be immobilised on the surface of dendrimer molecules. The dendrimer mole-
 40 cules are polymeric in nature, synthesised by the repetitive coupling of small building molecules. They are commercially available from Aldrich Chemicals in a range of molecular weights with a choice of terminal functional groups e.g. NH₂ or COOH. Heterobifunctional linkers can then be used in conjunction with conventional coupling chemistry to prepare the detecting label conjugates. For small hapten molecules, e.g. β-agonists, anabolic steroids or antibiotics, it is preferred that a small dendrimer (preferably no more than 16 surface groups) is coupled to a large dendrimer (typically more than 64 surface groups). The small molecular weight hapten (less than 1,000 Dalton) is coupled to the chemical groups on the small dendrimer followed by covalent attachment of the detecting label. The dendrimer conjugate may be purified by dialysis and gel permeation chromatography.

[0068] The test reagents contain multiple components (e.g. enzyme-labelled antibodies, fluorescent-labelled anti-
 bodies, latex-immobilised antibodies, dendrimer antibody-fluorophore conjugates, dendrimer antibody-fluorophore
 50 conjugates, dendrimer antibody-enzyme conjugates, enzyme-labelled haptens, fluorescent-labelled haptens, etc) as appropriate for particular panels of tests. The panels of tests possible are very diverse and can be chosen on the basis of clinical diagnosis (or veterinary diagnosis) as appropriate. For example, a desirable panel is for the detection of infectious diseases (e.g. hepatitis, HIV, syphilis, etc). Other panels include fertility hormones, cardiac markers, allergy proteins, etc. As well as clinical parameters, there is also the ability to detect large panels of drug residues.

[0069] The present invention permits the identification of individual compounds such as antibiotics. For example, a
 55 quantitative result can be obtained for up to 20 antibiotics on a device of surface area of 1 cm² simultaneously, in a time frame of minutes typically, with a sensitivity superior to that for HPLC/GCMS methods and comparable to that for conventional single parameter enzyme immunoassays. This approach may be easily extended to anabolic steroids,

beta-agonists, beta-blockers, pesticides, therapeutic drugs etc.

[0070] For analysis, chemiluminescence, bioluminescence or fluorescence may be suitable. The detection system is preferably a charge-coupled device (CCD) camera equipped to measure both fluorescent and chemiluminescent light. Briefly, the CCD camera collects the light signal generated from the test areas on the microfabricated device and converts this into relative light units (RLUs).

[0071] Fluorescent-based detection systems may be read directly, using appropriate optical filters for the labelling fluorophore.

[0072] A suitable chemiluminescent reagent is luminol, which can be analysed at a wavelength of 433-445 nm. Chemiluminescence may also be observed, based on detecting alkaline phosphatase-labelled biological molecules using 1,2-dioxetane.

[0073] To facilitate the detection of analytes, this invention preferably utilises a chemiluminescent detection system, using a CCD. A back-illuminated camera is preferred, to improve the capture efficiency at the wavelength of the light generated by the chemiluminescent light reaction (approximately 433-445 nm in the case of luminol).

[0074] The whole system may be operated by a personal computer where a specifically designed programme controls the X-Y table, dispenser unit, sample handling, temperature control, incubation times and the CCD camera. Figures 12-14 show the organisation of such a system.

[0075] Fig. 12 illustrates schematically the interaction of a personal computer (PC) having two control units 31,32. Unit 31 is in communication with a CCD imaging system represented at 33. Unit 32 is in communication with a dispenser unit and an X-Y translation table with sample tray represented at 34 and 35, respectively.

[0076] Fig. 13 is a schematic representation of the X-Y translation table. This drawing shows a sample platform 41 mounted on a linear actuator 42. X-Y translation is under the control of stepper motors 43,44 connected to respective drives 45,46. Translation is limited by "home position" sensors 47,48.

[0077] The sensitivity of labelled biological molecules and certain unlabelled biological molecules to light may make it necessary to perform the assays in the absence of light. The absence of light is achieved by constructing the case in a light-tight manner. The light-tight environment is also preferably temperature-controlled, e.g. within $\pm 0.2^\circ\text{C}$ or, preferably, $\pm 0.1^\circ\text{C}$, to ensure satisfactory assay precision and accuracy.

[0078] Fig. 14 shows the apparatus in perspective, part plan and cut-away side views. In particular, Fig. 14A shows reagent storage container 51, a light-tight door 52 and a camera body 53 with a removable cover 54. The major part of the camera can be external to the casing. The camera lens is placed in an aperture in the casing.

[0079] Fig. 14B shows in outline the samples on a sample tray 55 and, in outline, waste area 56 and an imaging area 57. The camera 53 is positioned over these areas. Fig. 14C shows, in addition to the container 51, camera 53, X-Y table 41 and stepper motor 43, a dispenser pump 58.

[0080] The design of the system shown in Fig. 14 is based on 3x3 sample rack holders of which 20 can be held at any one time. This means that, if 20 individual reaction regions are located on each 1 cm² of microfabricated device, a total of 3600 analyses may be performed simultaneously on a single sample. Alternatively, 180 samples may be analysed simultaneously for 20 different test parameters.

[0081] As indicated above, the analyte may be labelled. The ligand may also be labelled, allowing analysis by fractional occupancy.

[0082] The following Examples illustrate the invention.

Example 1 Sulphonamide Assay

[0083] In this example, 12 individual antibodies, each antibody specific for a single sulphonamide, were immobilised by covalent attachment by contact interactions on to discrete regions of a flat ceramic (aluminium oxide) substrate having a chemically modified surface. A multi-analyte assay was performed, using a competitive immunoassay format.

[0084] In more detail, ceramic substrates (1 cm x 1 cm) were ultrasonically cleaned using an alkaline detergent (RBS35, 5% v/v) followed by double deionised water and then placed in 6M HCl for 16 hours. The chips were then placed in chromic acid for 1 hour in an ultrasonic bath.

[0085] The substrates were washed exhaustively with double deionised water and acetone and then dried in an oven at 120°C for 2 hours. Following this pretreatment, the substrates were silanated using the organosilane γ -glycidoxypropyl trimethoxysilane (10% v/v) in anhydrous toluene, 4-dimethylaminopyridine (1.25 g/L) and triethylamine (1% v/v). This mixture was refluxed for 4 hours and then left overnight at room temperature. The substrates were washed with toluene and acetone before curing for 4 hours at 120°C.

[0086] Following the curing step, the substrates were placed in containers and stored at room temperature until required for spotting of sulphonamide antibodies. The sulphonamide antibodies were spotted using a BIODOT XY3000 dispenser. The 12 sulphonamides assayed were sulphadoxine, sulphamethizole, sulphachloropyridazine, sulphamethoxypropylpyridazine, sulphamerazine, sulphapyridine, sulphisoxazole, sulphathiazole, sulphamethazine, sulphaquinoxaline, sulphadimethoxine, and sulphadiazine.

[0087] Dispensed volumes of approx. 20 nl were employed for each sulphonamide antibody. The 12 sulphonamide antibodies which formed 12 discrete areas on the 1 cm² substrate were incubated for 2 hours at 37°C. The substrates were washed with phosphate-buffered saline (PBS) (pH 7.2) containing 2% casein (w/v) and then blocked in same buffer overnight at +2-8°C. After washing with PBS containing PEG300 (0.05% v/v), the devices were placed in a carrier.

[0088] Multi-sulphonamide standards (200 µl) and a cocktail of sulphonamide horseradish peroxidase conjugates (100 µl) were added to the reaction wells containing the device as appropriate and incubated for 15 minutes at room temperature. The standards contained 5, 10, 50 and 100 ng/ml for each of the 12 sulphonamides.

[0089] Thereafter the multi-sulphonamide devices were washed with PBS/PEG buffer to remove excess reagents and 300 µl chemiluminescent substrate [luminol (1.4mM)/urea hydrogen peroxide (9.6mM)] per device was introduced.

The devices were imaged using a CCD camera with an exposure time of up to 4 minutes. Standard curves were obtained for each of the 12 individual sulphonamides. Calibration curves for each of the 12 individual sulphonamides are represented graphically in Figure 15. The % B/B₀ value plotted on the Y-axis represents the % of the zero standard RLU (Relative light unit) value caused by each individual sulphonamide standard (plotted on the X-axis as log₁₀).

Example 2 Hormone Assay

[0090] In this example, a multi-analyte assay was performed for 3 hormones of large molecular weight, i.e. Prolactin (PL), Follicle-Stimulating Hormone (FSH) and Luteinising Hormone (LH). This example represents a multi-analyte assay for a sandwich-based immunoassay. No significant cross-reactivity was observed when the three hormones were determined in the same panel.

[0091] The chemical pretreatment and silanation procedures were exactly as described in Example 1. Individual PL, FSH or LH monoclonal antibodies (approx. 20 nl antibody dispensed) were immobilised on discrete areas of the chemically modified substrate. The multi-analyte assays were performed on both silicon and ceramic substrates with an epoxide surface as described in Example 1.

[0092] In the assay, 150 µl of a multiple LH/PL/FSH serum-based standard and 150 µl of a diluent assay buffer were added to the device and incubated for 15 minutes at room temperature. Following a wash step, 300 µl of a single conjugate cocktail of LH-HRPO/PL-HRPO/FSH-HRPO conjugates was added and incubated for 15 minutes. Thereafter the devices were washed to remove excess reagents and the chemiluminescent reagent [luminol (1.4 mM)/urea hydrogen peroxide (9.6 mM)] was introduced. The devices were imaged using a CCD camera with an exposure time of up to 4 minutes. Standard curves for each of the hormones were plotted after the images were processed.

Example 3 Sulphonamide Assay

[0093] In contrast to Example 1, a multi-sulphonamide assay has also been conducted using microchannels. The device is illustrated in Figure 11. In this Example, the reagent addition reservoirs 21 are 2 mm x 2 mm, and 300 µm deep (vol. 1.2 µl), the channels 23 are each 5 mm long, 200 µm wide and 100 µm deep (vol. 100 nl), and the reservoir 24 is 1.9 mm x 8.6 mm, and 300 µm deep (vol. 4.9 µl).

[0094] The chemical modification of the surface was performed as described in Example 1. Antibody was added to each of the channels and incubated for 2 hours at 37°C. The substrates were then blocked and washed as for Example 1.

[0095] A multi-sulphonamide standard (200 µl) and sulphonamide horseradish peroxidase conjugates (100 µl) were mixed. 1 µl of the resultant reagent was pipetted into each of the reservoirs supplying the antibody-coated channel for each sulphonamide. The standards contained 10 or 100 ng/ml of all the sulphonamides as appropriate.

[0096] The reagent flowed by capillary action. After incubation for 2 minutes, the substrate was washed 5 times with PBS/PEG and chemiluminescent reagent [luminol (1.4 mM)/urea hydrogen peroxide (9.6 mM)] was added.

[0097] The devices were imaged using a CCD camera with an exposure time of up to 4 minutes. The % B/B₀ values for the 4 sulphonamide curves are given in Table 1.

Table 1

Sulphonamide	%B/B ₀		
	0 ng/ml	10 ng/ml	100 ng/ml
Sulphamethazine	100	14	7
Sulphamethoxyypyridazine	100	28	15
Sulphaquinoxaline	100	56	24
Sulphamerazine	100	40	22

[0098] The utility of this invention by comparison with photolithography was demonstrated by the degree of non-specific adsorption of biological molecules on a photolabile (benzophenone-treated) substrate. The results are shown in Table 2.

5
10
15
Table 2

Mouse IgG	Irradiation By UV lamp (10 minutes)	Grey Mean (RLU)	
✓	✓	22368	24022
✓	X	17586	20531

20
25
[0099] Mouse IgG was detected using anti-mouse HRPO conjugate using chemiluminescent detection by CCD camera. The silicon or ceramic substrates having immobilised benzophenone photolabile linker should not bind mouse IgG when reacted in the absence of light. However, non-specific binding is occurring, since approximately 80% of the grey mean RLU achieved when the mouse IgG binding is performed under UV light is due to passive binding interactions. An array of biological molecules immobilised through covalent interactions according to this invention is demonstrably more distinct.

30
[0100] The evidence for covalent attachment is provided in the examples discussed below. In the first instance, ceramic substrates were silanated with APTES and then reacted with biotin-LC-sulpho NHS. The control ceramic substrate was not silanated with APTES; therefore, no terminal nucleophilic-NH₂ groups were available to react with the succinimidyl ester of the biotin derivative. The substrates were then reacted with an avidin-FITC conjugate and the fluorescence determined by CCD camera. Results are given in Table 3.

35
40
Table 3

[Biotin-LC-NHS]	APTES Substrate RLU (1 second exposure)	Control Substrate RLU (1 second exposure)
0	2,448	NS
100 µg/ml	8,881	NS
500 µg/ml	7,922	NS

40
NS=No fluorescent signal detected.

[0101] This clearly demonstrates specific immobilisation of biotin-LC-NHS. The maximum concentration of biotin-LC-NHS immobilised was approximately 100 µg/ml, since the RLU result for the 500 µg/ml substrate did not increase.

45
[0102] In a further example, ceramic substrates silanated with APTES were reacted with a dihydrazide linker. Sulphonamide antibodies were treated with sodium periodate to render them reactive towards the hydrazide linker. Control sulphonamide antibodies were simply dialysed against sodium acetate buffer pH 5.5. The RLU results are shown in Tables 4 (not treated) and 5 (treated by periodate method). The results clearly show that the hydrazide linker has been successfully linked to the ceramic surface. The control antibodies (no periodate activation) gave very poor standard curves when compared to the covalently immobilised sulphonamide antibodies.

50
[0103] The % of binding due to covalent or passive interactions is compared in Table 6. The RLU values for the zero and 10 ng/ml standards are provided in Table 7. Covalent interaction contributed on average 81.8% to the overall binding, clearly indicating covalent binding of the sulphonamide antibodies.

55
Table 4

SULPHONAMIDE	0 ng/ml	10 ng/ml	% B/B ₀	100 ng/ml	% B/B ₀
Sulphadoxine	2825	1456	51	NS	-

EP 0 874 242 B1

Table 4 (continued)

SULPHONAMIDE	0 ng/ml		10 ng/ml	% B/B ₀	100 ng/ml	% B/B ₀
Sulphamethizole	3531		1257	36	NS	-
Sulpachloropyridazine	6476		1585	24	NS	-
Sulphamethoxy-pyridazine	1099		928	84	NS	-
Sulphamerazine	2177		1137	52	1224	56
Sulphasoxazole	4879		1108	23	NS	-
Sulphathiazole	2932		814	28	NS	-
Sulphamethazine	NS		NS	-	NS	-
Sulphaquinoxaline	1041		828	79	968	93
Sulphadimethoxine	804		NS	-	781	97

Table 5

SULPHONAMIDE	0 ng/ml		10 ng/ml	% B/B ₀	100 ng/ml	% B/B ₀
Sulphadoxine	10520		3240	31	2135	20
Sulphamethizole	17141		5689	33	4882	28
Sulphachloropyridazine	24944		7565	30	2096	8
Sulphamethoxy-pyridazine	14082		10509	74	5687	40
Sulphamerazine	12594		5521	43	3240	26
Sulphasoxazole	24419		6686	27	2270	9
Sulphathiazole	14279		4602	32	2353	16
Sulphamethazine	3644		2810	77	2213	61
Sulphaquinoxaline	10575		6112	58	5588	53
Sulphadimethoxine	5526		2554	46	1983	36

Table 6

Sulphonamide	Percentage of Sulphonamide Antibody Binding	
	Covalent Interactions	Passive Interactions
Sulphadoxine	73.1	26.9
Sulphamethizole	79.4	20.6
Sulphachloropyridazine	74.0	26.0
Sulphamethoxy-pyridazine	92.2	7.8
Sulphamerazine	82.7	17.3
Sulphasoxazole	80.0	20.0
Sulphathiazole	79.4	20.6
Sulphamethazine	-	-
Sulphaquinoxaline	90.2	9.8
Sulphadimethoxine	85.4	14.6
Mean	81.8	18.2

EP 0 874 242 B1

Table 7

5	Sulphonamide	RLU			
		Covalent		Passive	
		0	10 ng/ml	0	10 ng/ml
	Sulphadoxine	32756	11131	1950	904
	Sulphamethizole	39020	11132	2782	1359
10	Sulphachloropyridazine	39632	8434	4410	1051
	Sulphamethoxy pyridazine	29489	13408	1793	770
	Sulphamerazine	28455	11077	2011	988
	Sulphisoxazole	38486	5774	4083	1031
	Sulphathiazole	28837	8087	2010	675
15	Sulphamethazine	11331	7535	802	574
	Sulphaquinoxaline	13838	8716	951	548
	Sulphadimethoxine	13062	5832	910	581
20	Covalent Immobilisation: Direct spotting on to glycidoxysilane surface.				
	Passive Immobilisation: Direct spotting on to dichlorodimethylsilane-reacted surface.				

[0104] Clearly, the covalent method results in superior results compared to the passive method. The results for the passive method can be increased approximately 2-fold by acid pretreatment of the sulphonamide antibodies, but this is still inferior to the covalent approach.

25 [0105] Confirmatory analysis on the chemically-modified silicon and ceramic substrates was also performed using X-ray photon spectroscopy (XPS). Survey spectra were recorded from two random areas of each substrate sample, from which their surface chemical compositions were determined. See Table 8 for results (given in atomic %).

Table 8

30	Sample	Area	C	O	Si	Al	N	Cl	Ca
	Silicon substrate (untreated)	1	21.9	52.7	25.4	-	-	-	-
		2	23.0	51.0	26.0	-	-	-	-
35	Silicon substrate silanated with APTES	1	55.5	23.3	10.5	-	10.2	0.5	-
		2	55.6	22.5	10.9	-	10.5	0.5	-
		1	51.3	25.6	13.0	-	7.2	2.9	-
		2	52.5	25.0	12.3	-	7.5	2.7	-
40	APTES silicon substrate treated with FITC	1	58.7	25.0	9.6	-	6.1	0.5	-
		2	58.8	24.7	9.8	-	6.0	0.8	-
45	Ceramic substrate (untreated)	1	27.2	46.3	11.3	13.3	-	1.4	0.5
		2	27.1	46.6	9.6	14.8	-	1.1	0.7
50	Ceramic substrate silanated with APTES APTES ceramic substrate treated with FITC	1	47.0	31.6	13.4	2.6	5.4	-	-
		2	45.9	31.7	13.7	3.3	5.3	-	-
		1	52.0	29.3	11.3	2.4	5.0	-	-
		2	51.0	30.6	11.7	2.3	4.5	-	-

[0106] The atomic composition results show a very good conversion of the parent silicon and ceramic substrates with APTES organosilane, with good reproducibility in surface composition indicated for the two areas tested on each sample. The FITC-labelled substrates showed 70% and 77% labelling of silicon and ceramic respectively.

55 [0107] Quantitative methods of fluorescence measurement on a dark silicon substrate have been compared with those on a white ceramic (aluminium oxide) substrate. The RLU results from CCD detection for fluorescent molecules (FITC) covalently linked to each substrate were compared to the quantitative method after the FITC molecules were stripped by the method of Hook *et al* (Langmuir 1991, Vol 7, 142-151).

Table 9

Substrate	CCD Exposure Time (sec)	Grey Mean		Number of Stripped FITC Molecules/Substrate
		Side 1	Side 2	
Ceramic	0.1	7483	7063	4.548×10^{15}
Silicon	10	753	612	1.412×10^{16}

[0108] The quantitative analysis of fluorescent label obtained by stripping FITC label from the substrate and measuring the signal by CCD camera (Table 9) shows that, despite the 1000-fold increase in signal of ceramic over silicon, there are actually significantly more FITC molecules present on the silicon substrate.

[0109] A further example of this phenomenon is shown by the results in Table 10, from a prolactin assay performed on silicon and ceramic substrates using fluorescent latex particles linked to a prolactin-detecting antibody as the detection system. RLU values are given (20 sec exposure).

Table 10

[Prolactin] Standard	RLU	
	Silicon Substrate	Ceramic Substrate
550 MIU/ml	814	8594
2200 MIU/ml	799	16735

[0110] The performance of ceramic provided superior results to silicon using this fluorescent detection method. Further, the problems of the dark body effect on silicon using fluorescence may be solved by employing chemiluminescence as the detection method. In a comparison between identical assays for FSH performed on silicon and ceramic, using chemiluminescent detection, the former required an exposure time 2-fold longer than that for ceramic to achieve the same RLU value.

[0111] In this specification, the following abbreviations apply:

APTES =	aminopropyltriethoxysilane
CK-MB =	creative kinase MB subunit
HRPO =	horseradish peroxidase
LC-sulfo-NHS =	long-chain sulfo-N-hydroxysuccinimide
FITC =	fluorescein isothiocarbamate

Claims

1. A method for the production of a solid state device for performing multi-analyte assays, comprising a substrate and a multiplicity of discrete reaction sites each bearing a ligand covalently bonded to the substrate, wherein the surface of the substrate between the reaction sites is inert with respect to analyte, the method comprising the steps of:

activating, and rendering hydrophobic, the said surface; and
applying to the activated, hydrophobic surface an array of ligands, to form the discrete reaction sites, wherein the ligands are applied in water.

2. A method according to claim 1, wherein the surface of the substrate is non-uniform, whereby an enhanced signal can be obtained from ligand-analyte interaction.

3. A method according to claim 2, wherein the surface of the substrate comprises an array of reaction channels, ridges, pillars, spots, chambers, dimples, wells or pits.

4. A method according to any preceding claim, wherein the substrate is of ceramic, glass, quartz or silicon.

5. A method according to any preceding claim, wherein the device has an area of less than 1 cm².

6. A method according to any preceding claim, wherein the area of each reaction site is less than 1 mm².
7. A method according to any preceding claim, which additionally comprises blocking the activated surface, between the reaction sites.
- 5 8. A method according to any preceding claim, wherein application of the ligands comprises an initial step of contacting the activated surface with an organosilane.
- 10 9. A method according to claim 8, wherein the organosilane has the formula (RO)₃Si-(CH₂)_n-X, wherein each R is a hydrocarbyl group, n is an integer, and X is a functional group.
11. A method according to claim 8 or claim 9, which includes the use of a bifunctional cross-linker to facilitate covalent attachment of biological ligands to the organosilane.
- 15 11. A method according to any of claims 8 to 10, wherein a photolabile cross-linker is used to react with the organosilane having a nucleophilic or electrophilic terminal group.
12. A method according to any of claims 1 to 7, which includes derivatisation of the surface with macromolecules such as polystyrene latex particles, dendrimers or polyethylene glycol containing chemical groups which facilitate covalent attachment of the ligands.
- 20 13. A method according to any preceding claim, which additionally comprises applying ligands that bind materials whose presence interferes with assaying an analyte.

25 **Patentansprüche**

1. Verfahren zur Herstellung einer Festphasenvorrichtung zur Durchführung von Multianalyt-Assays, die ein Substrat und mehrere diskrete Reaktionsstellen umfaßt, von denen jede einen kovalent an das Substrat gebundenen Liganden trägt, wobei die Substratoberfläche zwischen den Reaktionsstellen inert gegenüber dem Analyt ist, und das Verfahren folgende Schritte aufweist:
- 30 Aktivierung und Hydrophobierung der Oberfläche; und
Applikation einer Reihe von Liganden auf die aktivierte, hydrophobe Oberfläche, um die vereinzelt Reaktionsstellen zu bilden, wobei die Liganden in Wasser appliziert werden.
- 35 2. Verfahren nach Anspruch 1, in dem die Oberfläche des Substrates nicht gleichförmig ist, wodurch durch die Wechselwirkung zwischen Ligand und Analyt ein verstärktes Signal erreicht werden kann.
- 40 3. Verfahren nach Anspruch 2, in dem die Oberfläche des Substrates eine Anordnung von Reaktionskanälen, Gerten, Säulen, Flecken, Kammern, Erhebungen, Löchern oder Vertiefungen umfaßt.
4. Verfahren nach einem der vorstehenden Ansprüche, in dem das Substrat aus Keramik, Glas, Quarz oder Silizium ist.
- 45 5. Verfahren nach einem der vorstehenden Ansprüche, in dem die Vorrichtung eine Oberfläche von weniger als 1 cm² hat.
- 50 6. Verfahren nach einem der vorstehenden Ansprüche, in dem die Oberfläche jeder Reaktionsstelle weniger als 1 mm² beträgt.
7. Verfahren nach einem der vorstehenden Ansprüche, das außerdem eine Blockierung der aktivierten Oberfläche zwischen den Reaktionsstellen umfaßt.
- 55 8. Verfahren nach einem der vorstehenden Ansprüche, in dem die Applikation der Liganden einen einleitenden Schritt der Kontaktierung der aktivierten Oberfläche mit einem Organosilan umfaßt.
9. Verfahren nach Anspruch 8, in dem das Organosilan die Formel (RO)₃Si-(CH₂)_n-X hat, wobei jedes R eine Hydro-

EP 0 874 242 B1

carbylegruppe, n eine ganze Zahl, und X eine funktionellen Gruppe ist.

10. Verfahren nach Anspruch 8 oder 9, das die Verwendung eines bifunktionellen Crosslinkers umfaßt, um die kovalente Bindung biologischer Liganden an das Organosilan zu erleichtern.

5

11. Verfahren nach einem der Ansprüche 8 bis 10, in dem ein photolabiler Crosslinker verwendet wird, um mit dem Organosilan zu reagieren, das eine nucleophile oder elektrophile Endgruppe besitzt.

12. Verfahren nach einem der Ansprüche 1 bis 7, das eine Derivatisierung der Oberfläche mit Makromolekülen, wie etwa Polystyrol-Latex-Partikeln, Dendrimeren oder Polyethylen-Glykol-enthaltenden chemischen Gruppen umfaßt, welche die kovalente Bindung der Liganden erleichtern.

10

13. Verfahren nach einem der vorstehenden Ansprüche, das außerdem eine Applikation von Liganden umfaßt, die Materialien binden, deren Gegenwart den Assay eines Analyts beeinträchtigt.

15

Revendications

1. Procédé pour produire un dispositif en phase solide pour réaliser des déterminations sur plusieurs analytes, qui comprend un substrat et plusieurs sites de réaction séparés les uns des autres dont chacun comprend un ligand en liaison covalente avec le substrat, la surface du substrat étant inerte par rapport à l'analyte entre les sites de réaction et le procédé comprenant les étapes qui consistent à:

20

activer ladite surface et la rendre hydrophobe et

appliquer une série de ligands sur la surface activée et hydrophobe afin de former des sites de réaction séparés les uns des autres, les ligands étant appliqués dans de l'eau.

25

2. Procédé selon la revendication 1, dans lequel la surface du substrat n'est pas uniforme, ce qui permet d'obtenir un signal amélioré de l'interaction entre le ligand et l'analyte.

30

3. Procédé selon la revendication 2, dans lequel la surface du substrat comprend une série de canaux de réaction, d'arêtes de réaction, de colonnes de réaction, de taches de réaction, de chambres de réaction, de cuvettes de réaction, de puits de réaction ou de trous de réaction.

4. Procédé selon l'une quelconque des revendications précédentes, dans lequel le substrat est formé de céramique, de verre, de quartz ou de silicium.

35

5. Procédé selon l'une quelconque des revendications précédentes, dans lequel le dispositif présente une surface de moins de 1 cm².

40

6. Procédé selon l'une quelconque des revendications précédentes, dans lequel la surface de chaque site de réaction est inférieure à 1 mm².

7. Procédé selon l'une quelconque des revendications précédentes, qui comprend en outre l'étape qui consiste à bloquer la surface activée entre les sites de réaction.

45

8. Procédé selon l'une quelconque des revendications précédentes, dans lequel l'application des ligands comprend une étape initiale qui consiste à mettre la surface activée en contact avec un organosilane.

9. Procédé selon la revendication 8, dans lequel l'organosilane présente la formule (RO)₃Si-(CH₂)_n-X, dans laquelle chaque R représente un groupe hydrocarbyle, n représente un nombre entier et X représente un groupe fonctionnel.

50

10. Procédé selon la revendication 8 ou la revendication 9, qui comprend l'utilisation d'un agent de liaison bifonctionnel pour faciliter la liaison covalente de ligands biologiques avec l'organosilane.

55

11. Procédé selon l'une quelconque des revendications 8 à 10, dans lequel on fait réagir avec l'organosilane un agent photolabile de liaison qui comporte un groupe terminal nucléophile ou électrophile.

EP 0 874 242 B1

12. Procédé selon l'une quelconque des revendications 1 à 7, qui comprend la dérivatisation de la surface au moyen de macromolécules telles que des particules de latex de polystyrène, des dendrimères ou des groupes chimiques qui contiennent du polyéthylène glycol et qui facilitent la liaison covalente des ligands.

5 13. Procédé selon l'une quelconque des revendications précédentes, qui comprend en outre l'application de ligands qui se lient à des matériaux dont la présence interfère avec la détermination d'un analyte.

10

15

20

25

30

35

40

45

50

55

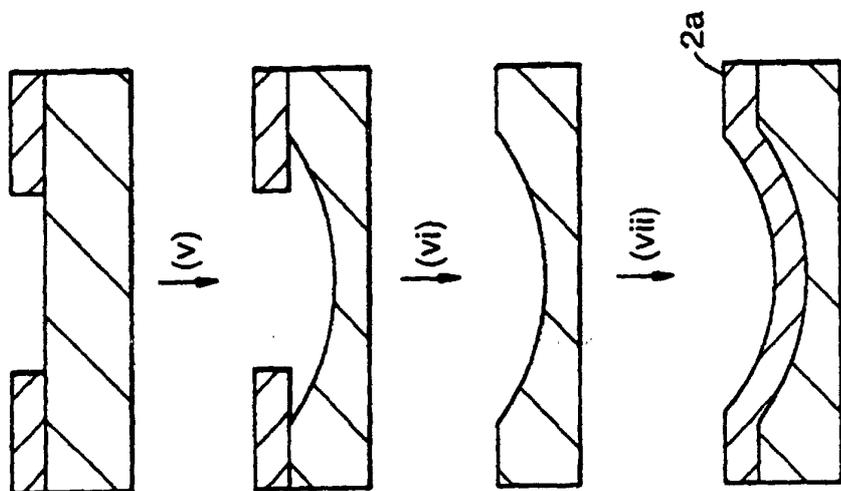


Fig. 1.

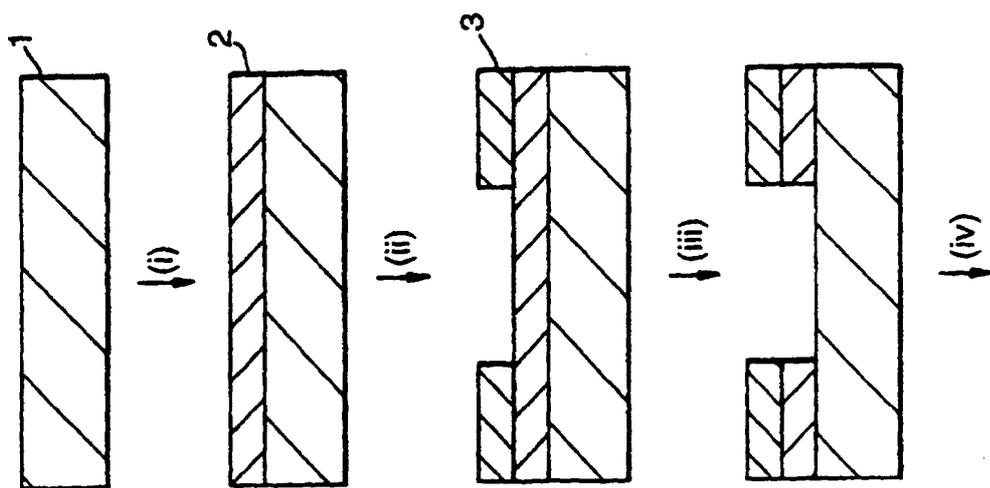


Fig.2.

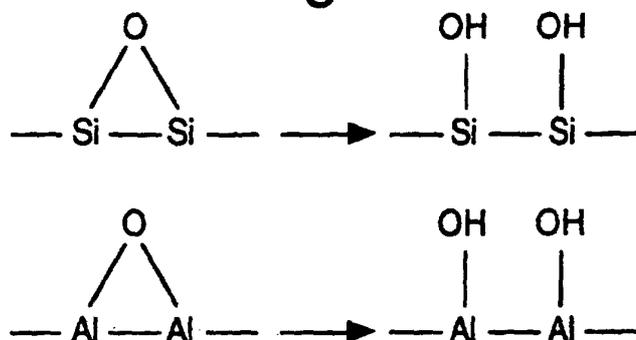


Fig.3.

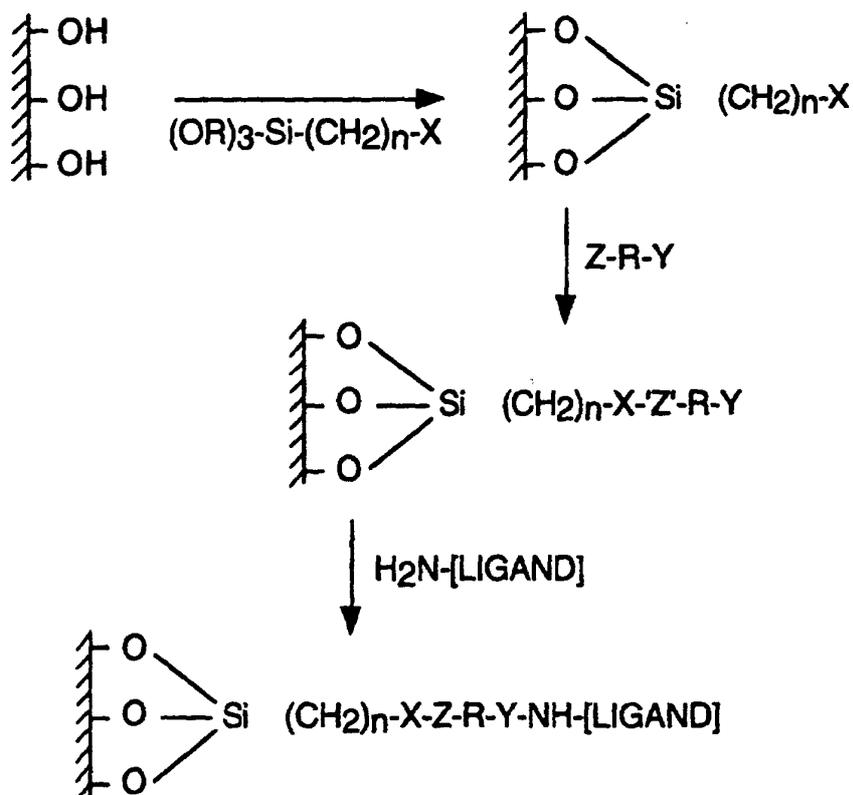


Fig.4.

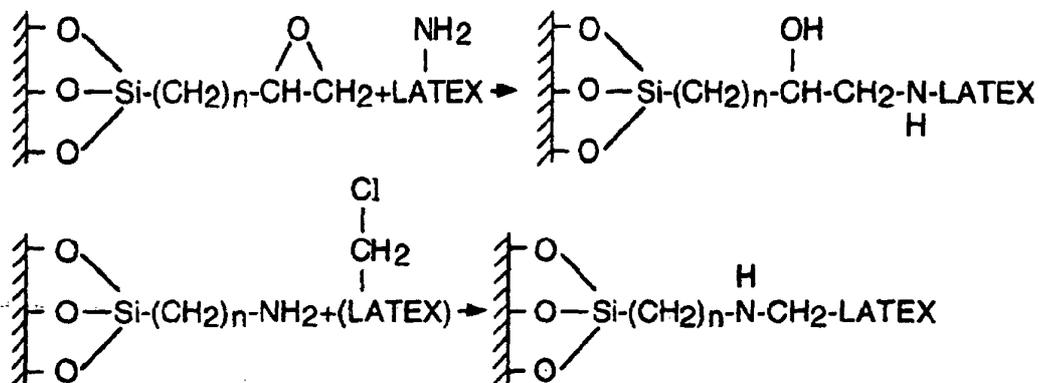


Fig.5.

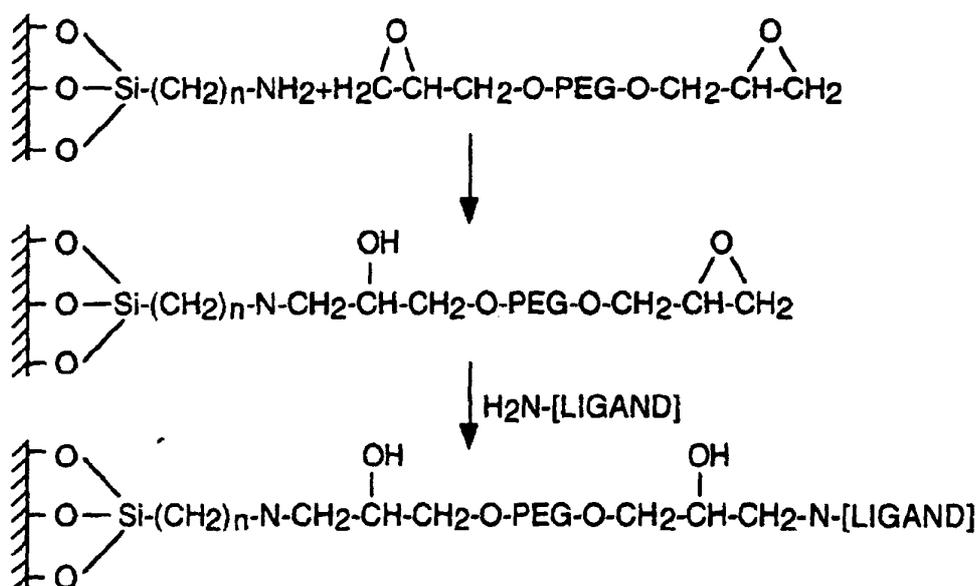


Fig.6.

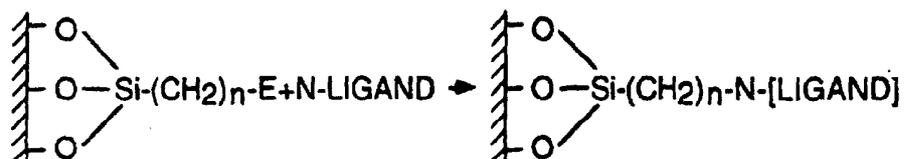


Fig.7.

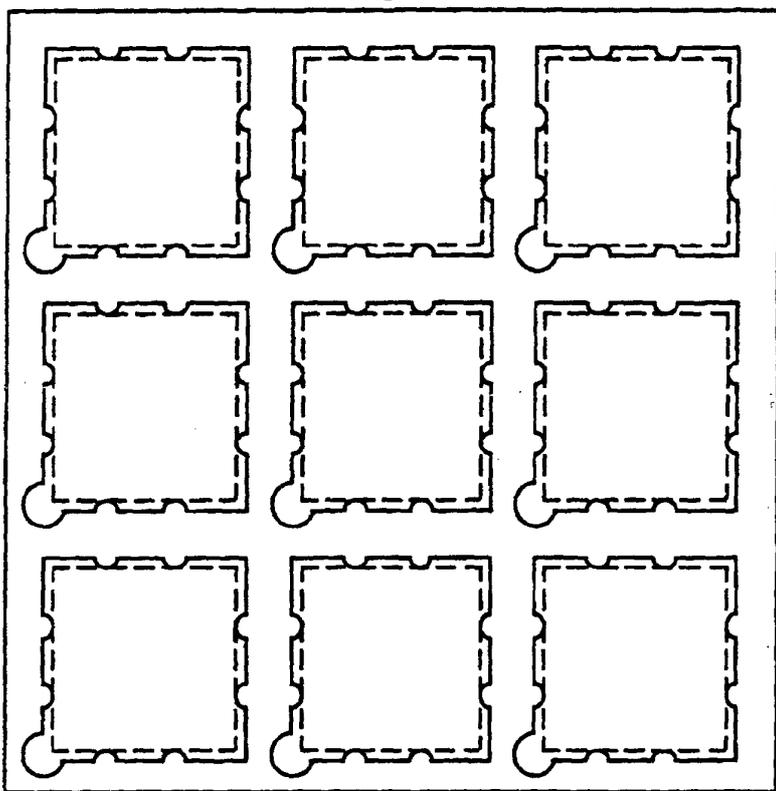


Fig.8.

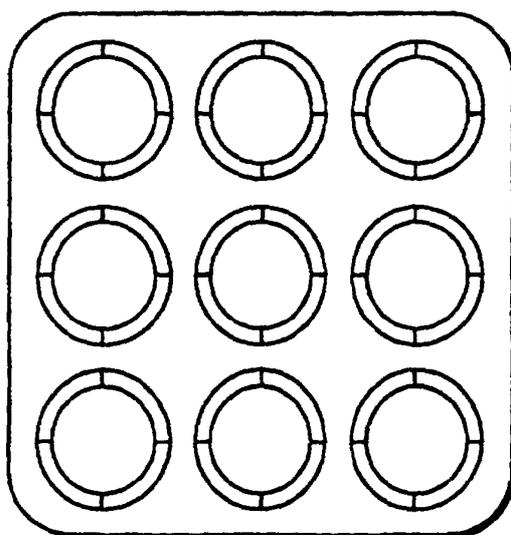


Fig.9.

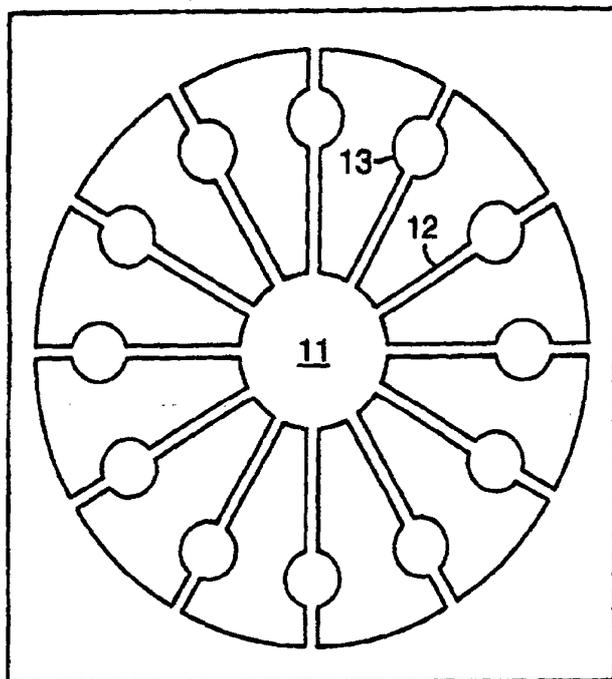


Fig.10.

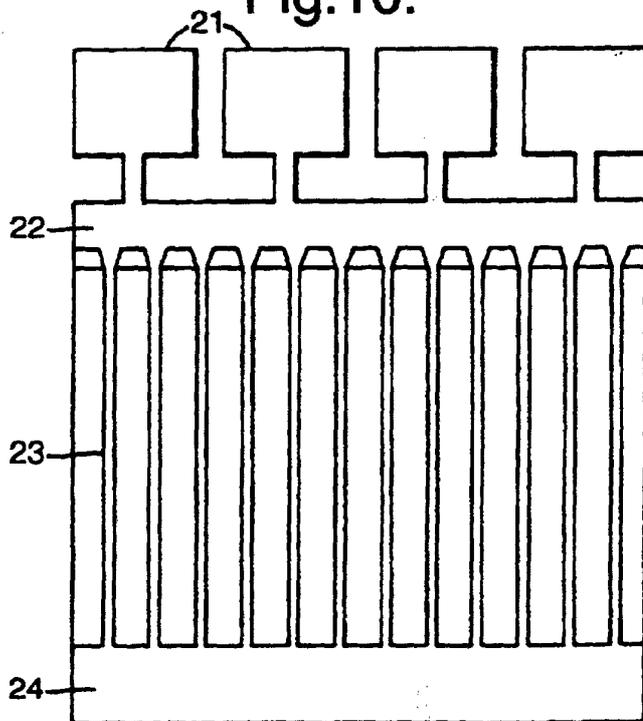
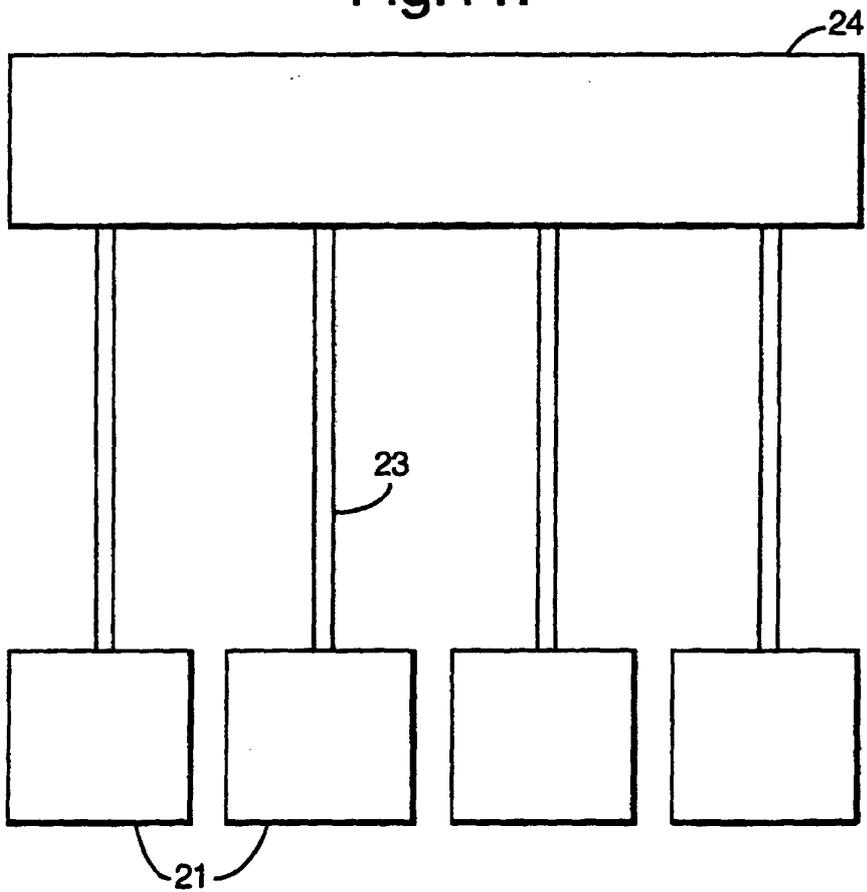


Fig.11.



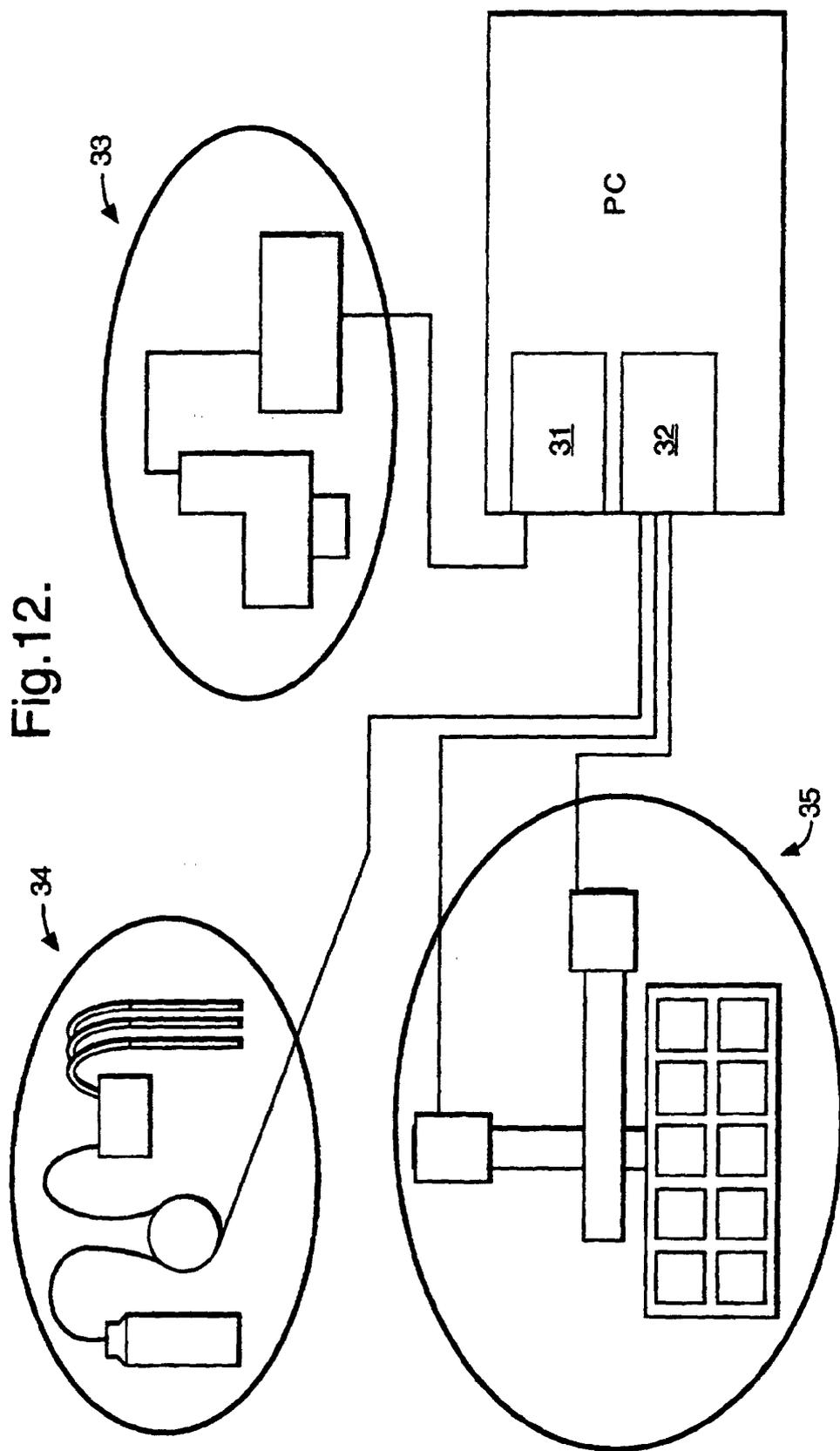


Fig.13.

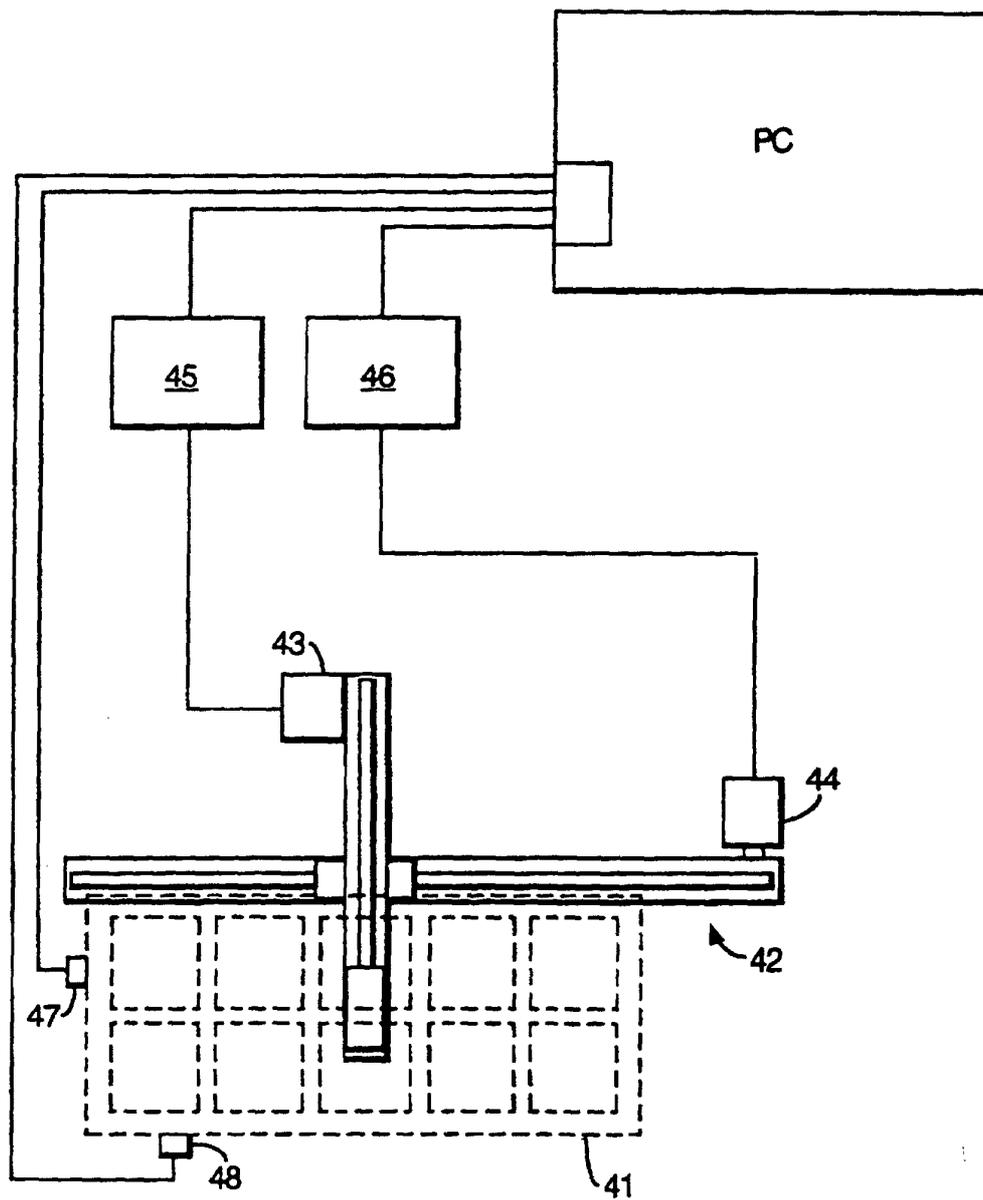


Fig.14A.

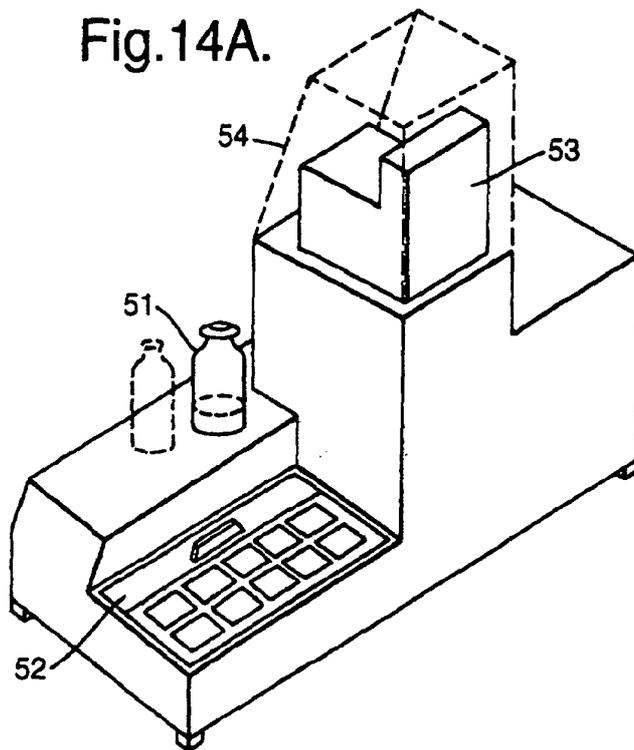


Fig.14B.

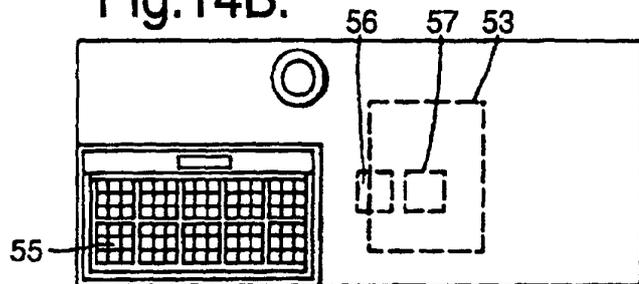


Fig.14C.

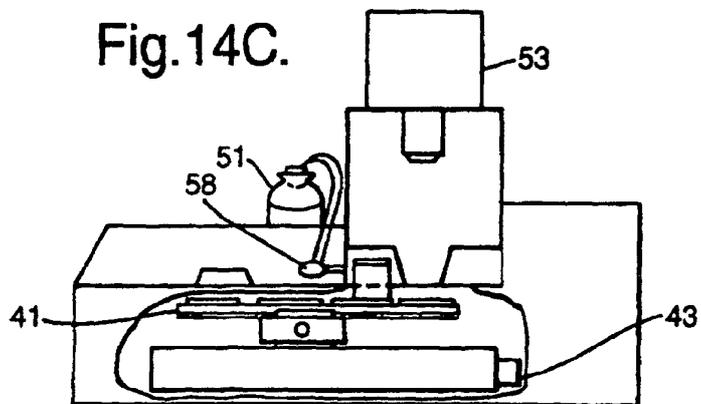


Fig.15.

